

App-Based Remote Alcohol Monitoring, Simplified



- ✓ Access video and location verified breathalyzer results from your smartphone
- ✓ Set random, recurring or custom testing schedules for individual users
- ✓ Receive real-time email or SMS notifications for failed and missed tests

COMPANY OVERVIEW

Founded: 2001

Number of Employees: 15

Headquarters: San Francisco, CA

Other Locations: San Diego, CA, Philadelphia, PA, Salt Lake City, UT

Founded in 2001, BACtrack is one of the pioneers of the breathalyzer industry. In 2004, it was the first company to obtain FDA approval to manufacture and market breathalyzers for consumer use in the United States. Since then, BACtrack has grown into a global leader in the alcohol testing and monitoring space, providing some of the most accurate and innovative breath alcohol testing devices and software on the market.


Launched in 2016, BACtrack View was designed to fit the specific needs of treatment centers, health clinics, businesses and individuals, making it easy and affordable to remotely monitor the sobriety of users, directly from a smartphone. It currently supports thousands of users, testing multiple times per day, typically several days per week. To date, over one million video breath alcohol tests have been recorded, processed and securely stored on the BACtrack View platform.

ABOUT THE SERVICE


BACtrack View is an app-based remote alcohol monitoring service, designed to track the sobriety of users by providing a video of their breath alcohol tests to the person or organization monitoring them.

Testers are given a BACtrack Mobile Smartphone Breathalyzer to use while on the service, which connects to the BACtrack View App on an iPhone or iPad via Bluetooth. When it's time to take a test, the Tester receives a notification on their phone, which prompts them to turn on their breathalyzer and record themselves taking the test. The BAC results and video, along with the date, time and GPS location are uploaded to the BACtrack View App in near real-time for review by their Manager

Each Tester is assigned a Manager, who can request on-demand tests, set recurring or randomized testing schedules, view previous test video results, and download a summary PDF report of that user's test history. Managers have the option to receive Email or SMS notifications when an MU fails or misses a test. Additional Managers within the organization also have the option to receive these notifications.

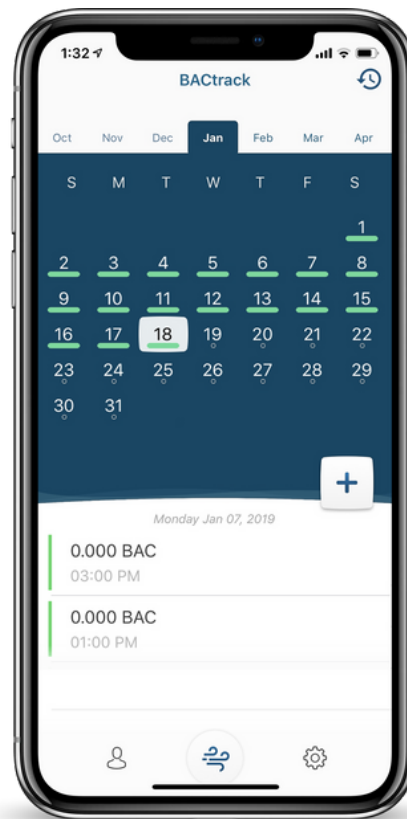

HOW IT WORKS

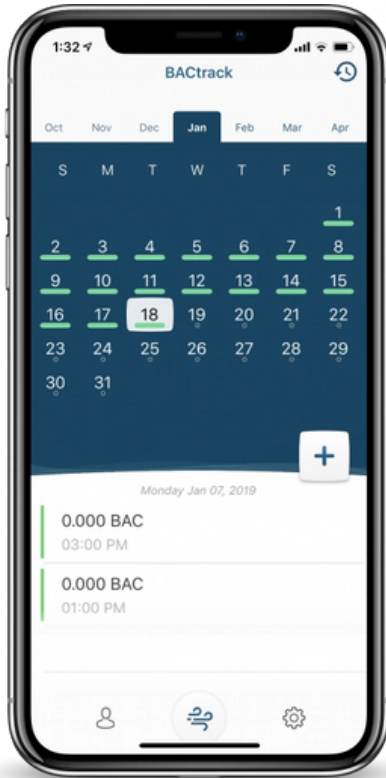


Set Testing Schedule

Choose on-demand or schedule random alcohol testing for someone you'd like to monitor remotely or for self-testing.



Request a BAC Test

Automated notifications tell a user when it's time to test. BACtrack Mobile and the BACtrack View app offer easy and discreet BAC testing.

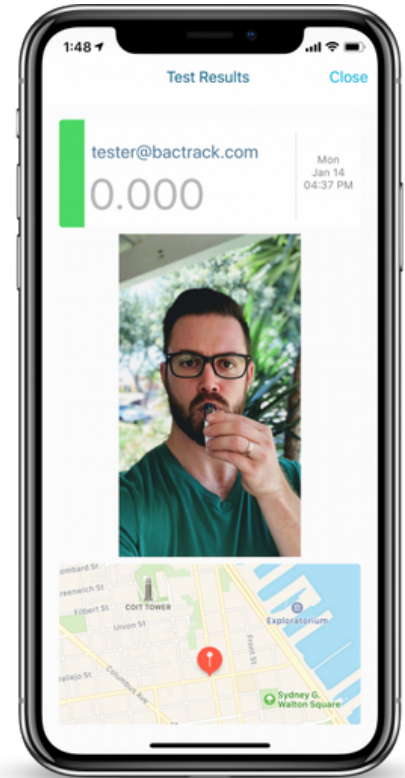

Get Real-Time Results

Video and location-verified BAC results are captured and saved via the App, and automatically shared with Monitors via Text or Email.

SERVICE SUMMARY

- **Smartphone App for Managers**

- Request on-demand BAC tests
- Set recurring or randomized testing schedules
- View video BAC results
- Receive real-time test result notifications via text or email
- Download PDF testing history

- **Smartphone App for Monitored Users**

- Text notifications when it's time to take a test
- Records video of users during their test
- Can also take on-demand BAC tests at any point

- **High quality, discrete, fuel-cell breath alcohol tester**

- Bluetooth connectivity to App
- Charges via micro USB cable (included) - 200-300 tests on one battery charge
- 2 washable and reusable mouthpieces included, required for proper test results

CURRENT & PAST PARTNERS

- Washoe County, NV
- Warrick County, IN
- Bayfield County, WI
- Denton County, TX
- Lane County, OR
- Brown County, WI
- Mahoning County, OH
- Clayton County, GA
- Benton County, AR
- Kent County, MI
- Killeen County, TX
- San Diego County, CA
- Atascosa County, TX
- Calhoun County, MI
- Bell County, TX
- Clark County, NV
- Fresno County, CA
- Berkeley County, SC
- Davis County, UT
- Dekalb County, GA
- Grants Pass, OR
- Orleans County, NY
- Richland County, SC
- St. Louis County, MN
- Warrick County, IN
- Waterford, MI
- Woodbury, MN
- Tiverton, RI
- Rochester, MN
- Newcastle, DE
- Lynchburg, VA
- Manchester, CT
- Alexandria, VA
- Bartholomew, IN
- Bayfield County, WI
- Brazos County, TX
- Cass County, ND
- Pleasanton, CA
- Chochise, AZ
- Coon Rapids, MN
- Cook County, IL

PRICING

Basic Apple iOS Only	Plus Apple iOS Only Popular	Pro Apple iOS Only
<p style="text-align: center;">\$ 79.99 / month</p>	<p style="text-align: center;">\$ 99.99 / month</p>	<p style="text-align: center;">\$ 129.99 / month</p>
<p style="text-align: center;">FREE* Pro-Grade Smartphone Breathalyzer [?]</p>	<p style="text-align: center;">FREE* Pro-Grade Smartphone Breathalyzer [?]</p>	<p style="text-align: center;">FREE* Pro-Grade Smartphone Breathalyzer [?]</p>
<p style="text-align: center;">30 Tests Per Month [?]</p>	<p style="text-align: center;">Unlimited Tests Per Month [?]</p>	<p style="text-align: center;">Unlimited Tests Per Month [?]</p>
<p style="text-align: center;">None PDF Reports [?]</p>	<p style="text-align: center;">1 Per Month PDF Reports [?]</p>	<p style="text-align: center;">Unlimited PDF Reports [?]</p>
<p style="text-align: center;">1 Connections Allowed [?]</p>	<p style="text-align: center;">3 Connections Allowed [?]</p>	<p style="text-align: center;">Unlimited Connections Allowed [?]</p>
<p style="text-align: center;">Text Monitor Notifications [?]</p>	<p style="text-align: center;">Text Monitor Notifications [?]</p>	<p style="text-align: center;">Text & Email Monitor Notifications [?]</p>
<p style="text-align: center;">Email Customer Service Level</p>	<p style="text-align: center;">Phone & Email Customer Service Level</p>	<p style="text-align: center;">Dedicated Account Rep Customer Service Level</p>
<div style="text-align: center; border: 1px solid #ccc; padding: 10px; margin: 10px auto; width: 80%;"> <p>Start 14 Day Free Trial</p> <p style="font-size: 0.8em;">Risk Free, Cancel Anytime!</p> </div>	<div style="text-align: center; border: 1px solid #0056b3; padding: 10px; margin: 10px auto; width: 80%; background-color: #0056b3; color: white;"> <p>Start 14 Day Free Trial</p> <p style="font-size: 0.8em;">Risk Free, Cancel Anytime!</p> </div>	<div style="text-align: center; border: 1px solid #ccc; padding: 10px; margin: 10px auto; width: 80%;"> <p>Start 14 Day Free Trial</p> <p style="font-size: 0.8em;">Risk Free, Cancel Anytime!</p> </div>